



Coding and Robotics
Applied Research Department

SMART Platform Imaging Software Operations Manual

Contents

Introduction	3
1 Hardware and Software	4
1.1 Hardware	4
1.2 Software	5
1.2.1 System Requirements	5
1.2.2 Installation	6
1.2.3 Table Operation	7
1.2.4 Data Processing	10
2 Output Files	16
3 Experiment Setup Recommendations	17
3.1 Background	17
3.2 Substrate and Container Visibility	18
3.3 Lighting	19
3.4 Plant Suitability	20
3.5 Pre-Experiment Setup Checklist	21
4 Troubleshooting	22
5 Future Considerations	23
5.1 Hardware	23
5.2 Software	23

Introduction

THIS manual describes the installation and operation of the SMART Platform Imaging Software. The software is designed to process experimental imaging datasets generated by the SMART platform, and to consolidate the resulting data into structured Excel reports for visualization and group analysis.

The software consists of two primary modules:

- **Analysis** — Processes raw image folders into motion and size data, outputting per-plant Excel spreadsheets.
- **Data Consolidation** — Merges all per-plant Excel files into a single workbook, with group assignment, chart generation, and optional plant filtering.

Before installation or operation, verify that the target computer satisfies the minimum hardware requirements outlined in Chapter 1.

Chapter 1

Hardware and Software

This chapter describes the physical SMART platform hardware and the software used to process the imaging data it produces.

1.1 Hardware

The SMART platform consists of a motorized X/Y gantry that positions a camera over each plant in turn for imaging. The primary hardware components are listed below:

- **Three NEMA 17 stepper motors** — Three 1.8° per step motors provide motion along the X, Y, and Z axes.
- **GRBL-based CNC controller** — Controls all motion by interpreting G-code commands received over a serial connection.
- **Limit switches** — Installed at the machine home position to establish the origin during the homing procedure.
- **24 V, 14.6 A power supply** — Provides power to the motion control electronics and stepper motors.
- **USB camera** — Mounted to the gantry for image acquisition at each plant location.
- **Light Strips** — Approximately a minimum of 15 square feet of lighting laid across the top of the table to ensure lighting quality standard for imaging.

The GRBL controller communicates with the SMART Table Settings software over a serial connection at 115200 baud. During operation, the machine is

first homed using the limit switches before executing the programmed imaging sequence. The camera captures an image at each predefined position as the gantry moves across the table.

Note: The limit switches are located only at the home (origin) position. No limit switches are installed at the maximum travel extents, so movements beyond the machine's physical limits are not automatically prevented once the machine has left the home position.

1.2 Software

This section covers the SMART Platform Imaging Software: its system requirements, installation procedure, operation of the imaging table, and the data processing pipeline.

1.2.1 System Requirements

Before installing the SMART Platform Imaging Software, ensure that the computer meets the minimum hardware specifications listed below.


Requirement	Minimum Specification
Processor	8 CPU cores or greater
Available Memory	12 GB RAM or greater
Free Storage Space	10 GB minimum, or twice the size of the dataset being processed
Operating System	Windows 10 or Windows 11

The processor and memory specifications can be verified using the Windows *System Information* utility.

Checking System Specifications

1. Press the **Windows** key.
2. Type *System Information*.
3. Open the *System Information* application.
4. Select *System Summary* from the navigation pane on the left.

The navigation pane should appear similar to the example shown below.

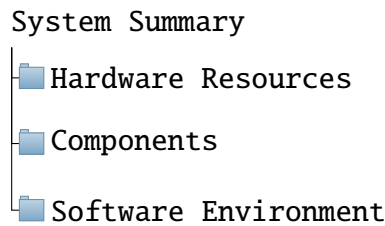

The *System Summary* page contains information similar to the table below.

Item	Value
Processor	...
Installed Physical Memory (RAM)	XX GB
Total Physical Memory	XX GB
Available Physical Memory	XX GB

Verify that:

- The processor contains at least **8 CPU cores**.
- The **Available Physical Memory** value is at least **12 GB**.
- The storage device used for processing has sufficient free space.

The recommended free storage capacity is either:

- At least **10 GB**, or
- At least **twice the size** of the experimental dataset being processed,

whichever value is greater.

WARNING: If the minimum hardware requirements are not met, the software may become unresponsive or terminate unexpectedly during processing. This may result in loss of processed data and, in some cases, corruption of experimental datasets.

1.2.2 Installation

This section describes the installation procedure for the SMART Platform Imaging Software.


To download the software onto the target machine, proceed to the [Github Repository](#). Once on this page, download both the **mergeData** and **Analysis** folders. The installation directory is arbitrary, as long as the user has read and write permissions to the chosen location.

Once the folders have been downloaded, no further installation steps are required. The software can be launched directly from its folder.

You have now installed the SMART Platform Imaging Software onto your machine.

1.2.3 Table Operation

Physical operation of an imaging run — moving the gantry over the table and capturing photos — is controlled through a separate **SMART Table Settings** interface. This interface connects to the table's motion controller and camera, drives the gantry through a fixed imaging pattern, and captures a photo at each plant position.

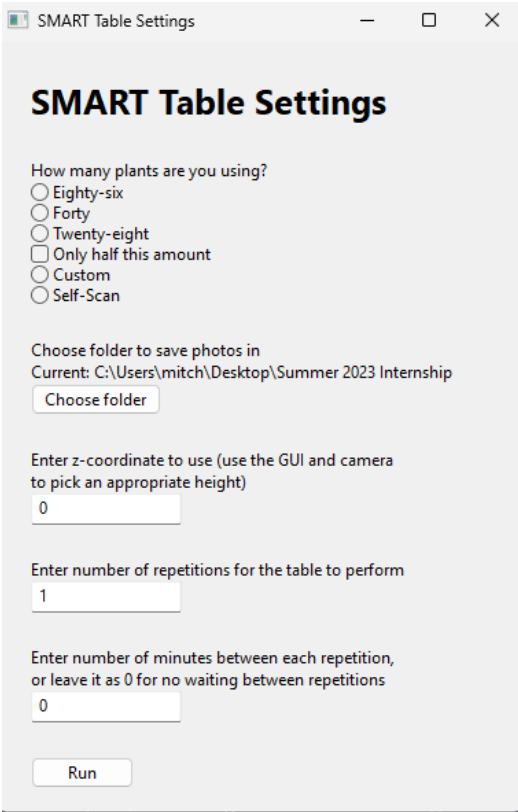

The screenshot shows a window titled "SMART Table Settings" with a standard Windows title bar (minimize, maximize, close buttons). The window content includes a title "SMART Table Settings" in bold. Below the title, there is a section "How many plants are you using?" with five radio button options: "Eighty-six", "Forty", "Twenty-eight", "Only half this amount", "Custom", and "Self-Scan". The "Only half this amount" option is selected. Below this is a section "Choose folder to save photos in" with the text "Current: C:\Users\mitch\Desktop\Summer 2023 Internship" and a "Choose folder" button. Below that is a section "Enter z-coordinate to use (use the GUI and camera to pick an appropriate height)" with a text input field containing "0". Below that is a section "Enter number of repetitions for the table to perform" with a text input field containing "1". Below that is a section "Enter number of minutes between each repetition, or leave it as 0 for no waiting between repetitions" with a text input field containing "0". At the bottom of the window is a "Run" button.

Figure 1.1: SMART Platform Control Software

Prerequisites

Before running the table for the first time, a `config.txt` file must exist in the program folder containing three lines, in order:

1. The destination folder where photos will be saved.
2. The serial port name used by the table's controller (e.g. COM3).
3. The numeric index of the USB camera to use.

The destination folder can be changed at any time from within the interface (see below); the serial port and camera index must currently be set by editing `config.txt` directly.

Launching the Interface

Run the SMART Table Settings program. A window titled **SMART Table Settings** will open, containing the controls described below.

Selecting a Layout

Choose the imaging pattern to use, based on how the plants are arranged on the table:

- **Eighty-six** — 86-plant grid layout.
- **Forty** — 40-plant grid layout.
- **Twenty-eight** — 28-plant grid layout.
- **Custom** — Uses a fixed list of coordinates read from `config.txt` instead of a built-in grid pattern.
- **Self-Scan** — Automatically locates plants by detecting green pixels through the camera, then records their coordinates to `config.txt` for later use with Custom mode.

Check **Only half this amount** to image only half of the positions in the selected grid layout (43, 20, or 14 plants respectively). This option has no effect in Custom or Self-Scan mode.

Choosing a Photo Folder

Click **Choose folder** to select the destination folder for captured images. The currently selected folder is shown above the button and is saved to `config.txt` automatically.

Setting Camera Height

Enter a Z-coordinate in the height field to set the camera height used for every photo in the run. Determine an appropriate value beforehand using the table's positioning controls together with the camera preview. This is done in the `config.txt` file.

Repetitions and Interval

- **Repetitions** — The number of times the table should image every plant, i.e. how many photos of each plant will be captured over the course of the run.
- **Interval** — The number of minutes to leave between the start of each repetition. Leave at **0** to run repetitions back-to-back with no wait. If a repetition takes longer than the specified interval to complete, the next repetition begins immediately.

Running the Table

Click **Run** to begin. The interface will:

1. Connect to the table's motion controller over serial and initialize the camera.
2. Home the machine.
3. Move to each plant position in turn, capturing and saving a photo at each stop.
4. Repeat the process for the configured number of repetitions, waiting between repetitions if an interval was set.

Once **Run** is clicked, the settings window closes and progress messages are printed to the console, including an estimate of the total run time when an interval is set.

WARNING: Do not disconnect the camera or the controller, or interrupt power to the table, while a run is in progress. Doing so may leave the gantry in an unknown position and require re-homing before the table can be used again.

NOTE: The table's limit switches are located only at the home (origin) position. Movement commands that would drive the gantry out of bounds are not automatically caught away from home, so custom coordinates and camera heights should be verified carefully before an unattended run.

1.2.4 Data Processing

This section describes the data processing side of the SMART Platform Imaging Software. The software has two primary functions:

1. **Analysis** — Processes raw experimental image datasets and generates per-plant Excel reports containing motion, size, and supplementary data.
2. **Data Consolidation** — Merges the per-plant Excel reports produced by the Analysis function into a single unified workbook, with support for group assignment, chart generation, and selective plant filtering.

Both functions are accessible from within the same application window via separate tabs at the top of the interface.

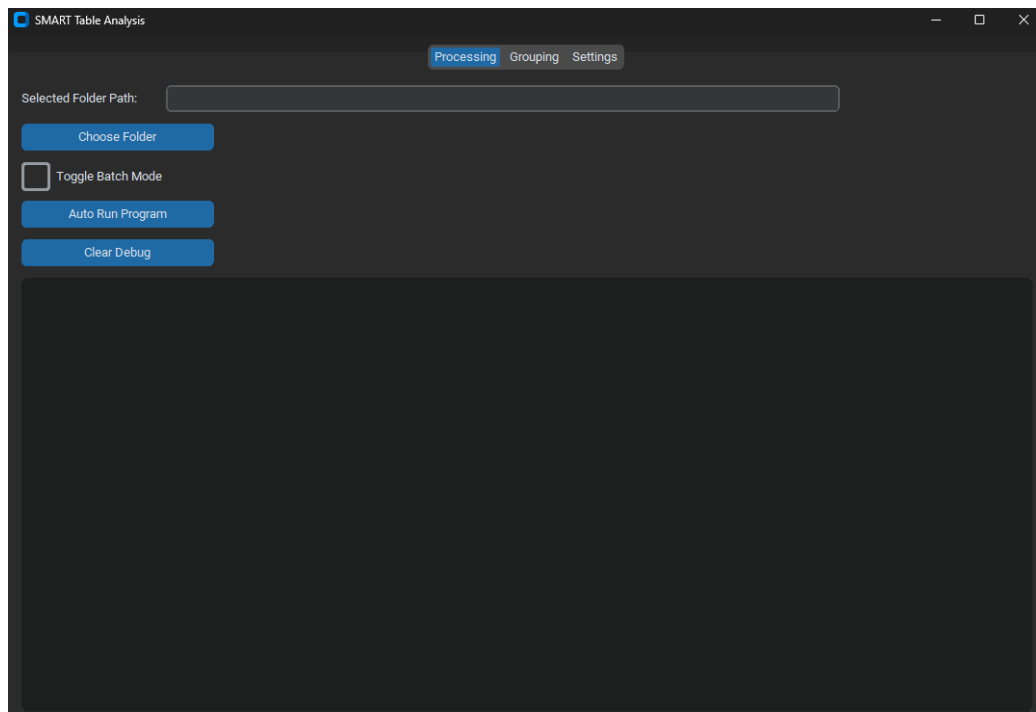


Figure 1.2: Processing Tab

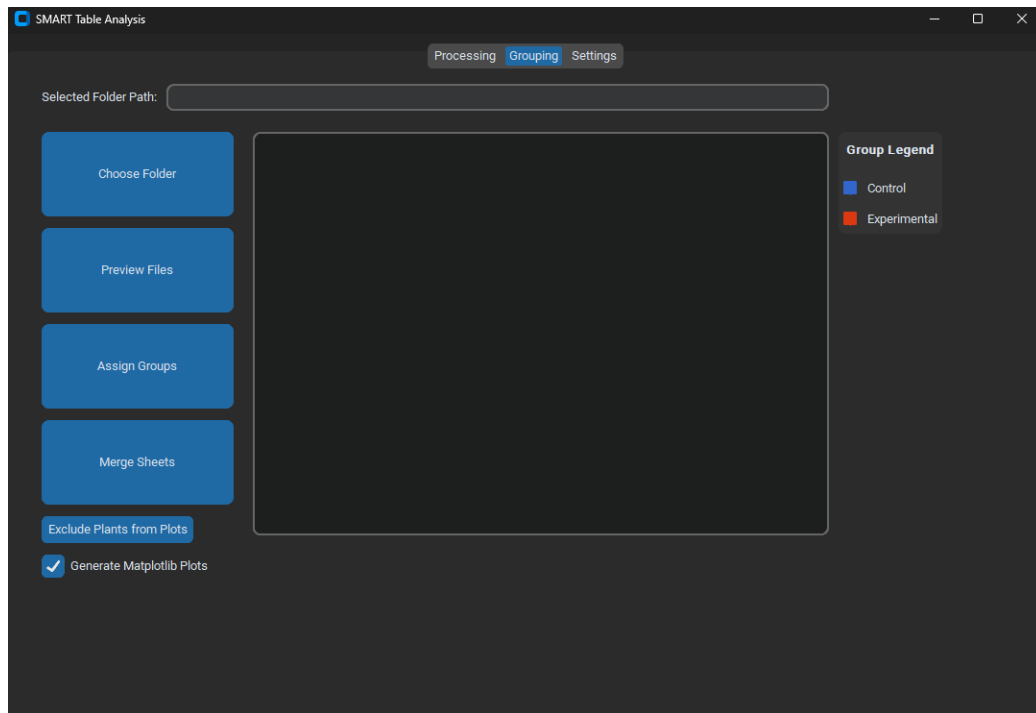


Figure 1.3: Grouping Tab

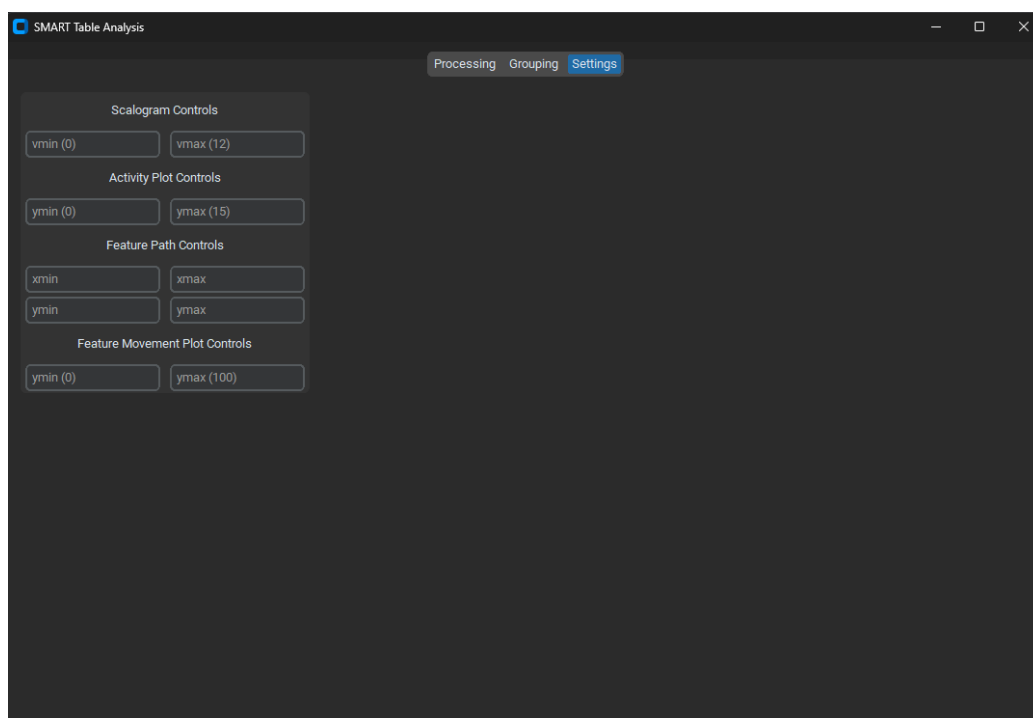


Figure 1.4: Settings Tab


Analysis

The Analysis function processes datasets and generates an Excel spreadsheet containing movement, size, and other derived data for each plant. This spreadsheet will be saved within the folder that was selected for processing.

Processing Modes There are two processing modes available:

1. **Single Folder Mode** — Used to process a single plant folder, for example Plant1. Select the folder containing the raw image files directly.
2. **Batch Processing Mode** — Used to process an entire experiment at once. The selected folder should contain one subfolder per plant, following the topology shown below.

SL_Tom_11-8-24

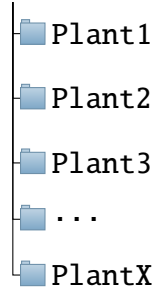

Selecting a Processing Mode To select a processing mode, check or uncheck the **Toggle Batch Mode** checkbox in the Processing tab. When the checkbox is unchecked, the program operates in Single Folder Mode. When checked, the program operates in Batch Processing Mode. The button label updates to reflect the current mode.

Running the Analysis To run the Analysis function, follow the steps below:

1. Select the desired processing mode using the **Toggle Batch Mode** checkbox.
2. Click the **Choose File** or **Choose Folder** button and navigate to the target folder.
3. Click **Auto Run Program** to begin processing.

Progress and status messages will appear in the debug output window below the buttons. If an error occurs during processing, an error message will be displayed in the debug window along with the name of the file or folder that caused the issue.

WARNING: Do not close the application or navigate away from the Processing tab while the program is running. Doing so may interrupt processing and result in incomplete or missing output files.

When processing is complete, the output Excel file will be saved inside the processed plant folder. For batch mode, one Excel file will be generated per plant subfolder.

Data Consolidation

The Data Consolidation function is accessed from the **Grouping** tab. It collects all Excel files produced by the Analysis function from a selected folder, merges them into a single workbook, assigns plants to experimental groups, and generates charts for visualization.

Step 1: Select a Folder Click the **Choose Folder** button and navigate to the root directory of the experiment. This should be the same parent folder used during batch processing, containing all plant subfolders.

Step 2: Preview Files Click **Preview Files** to scan the selected folder for Excel files. The debug window will display a list of all detected files and the number of plants identified. This step must be completed before assigning groups or merging data.

Step 3: Assign Groups Click **Assign Groups** to open the group assignment window. Each detected plant will be listed alongside radio buttons corresponding to the available experimental groups.

- By default, two groups are available: **Control** and **Experimental**.
- To add a new group, type a name into the text field at the bottom of the window and click **Add Group**.
- Select the appropriate group for each plant using the radio buttons.
- Click **Done** when all assignments are complete.

The group legend on the right side of the Grouping tab will update to reflect all current groups and their assigned colors.

Step 4: Exclude Plants from Plots (Optional) Click **Exclude Plants from Plots** to open the exclusion window. Plants are listed with checkboxes. Uncheck any plant to exclude it from the generated Matplotlib charts. Plants excluded here will still be included in the Excel output.

Step 5: Configure Chart Generation (Optional) Check or uncheck the **Generate Matplotlib Plots** checkbox to enable or disable the generation of .png chart files alongside the Excel output. This option is enabled by default.

Step 6: Merge Sheets Click **Merge Sheets** to begin the consolidation process. The program will:

1. Read all detected motion and size Excel files.
2. Merge the data into a single workbook with two worksheets: **Motion_Worksheet** and **Size_Worksheet**.

3. Generate embedded Excel charts for both motion and size data, color-coded by group.
4. Optionally generate external .png chart files (motion_chart.png and size_chart.png).
5. Save the output as Merged.xlsx in the selected folder.

A confirmation dialog will appear when the file has been created successfully. If Merged.xlsx already exists in the selected folder, the program will prompt you to confirm before overwriting it.

WARNING: Do not modify or open any of the source Excel files while the merge is in progress. This may cause file access errors and prevent the output from being generated.

Settings

The **Settings** tab provides controls for adjusting the axis limits of the plots generated during the Analysis phase. The available controls are listed below.

Control	Description
Scalogram vmin / vmax	Minimum and maximum color scale values for scalogram plots
Activity Plot ymin / ymax	Y-axis bounds for activity line plots
Feature Path xmin / xmax	X-axis bounds for feature path plots
Feature Path ymin / ymax	Y-axis bounds for feature path plots
Feature Movement ymin / ymax	Y-axis bounds for feature movement plots

Leave any field blank to use the default value shown in the placeholder text. Settings are applied the next time **Auto Run Program** is executed.

Chapter 2

Output Files

This chapter describes the files generated by the software and their locations.

Analysis Output

For each processed plant folder, the following files are generated inside that folder:

File	Description
<code>motion_analysis_report.xlsx</code>	Contains per-frame displacement data for the plant on the <code>Displacement Data</code> sheet.
<code>mask_sizes.xlsx</code>	Contains per-frame pixel count data for the plant on <code>Sheet1</code> .

Data Consolidation Output

The following files are generated inside the root experiment folder after running **Merge Sheets**:

File	Description
<code>Merged.xlsx</code>	Consolidated workbook containing motion and size data for all plants, with embedded Excel charts color-coded by group.
<code>motion_chart.png</code>	Matplotlib chart of plant displacement over time (generated if the plot option is enabled).
<code>size_chart.png</code>	Matplotlib chart of plant pixel count over time (generated if the plot option is enabled).

Chapter 3

Experiment Setup Recommendations

THE SMART Platform Imaging Software relies on color-based image segmentation to isolate plants from their surroundings. The accuracy of this segmentation — and therefore the reliability of all motion and size data — is directly determined by the physical setup of the experiment. This chapter describes the conditions under which segmentation performs best, and the setup choices that should be avoided.

The Platform should be in a air controlled environment. The electronics will give off heat sufficient to raise the ambient air temperature depending on the size of the room.

NOTE: The segmentation algorithm is designed for **majority-green plants**. Species or growth stages where the plant is predominantly yellow, brown, purple, or white are not reliably segmented and may produce inaccurate size measurements.

3.1 Background

The single most important factor in segmentation quality is the contrast between the plant and its background. The software identifies plant pixels by detecting green hues that stand out from the surrounding area. If the background contains similar colors or textures, the segmentation boundary becomes ambiguous.

Recommended

- Use a **uniform dark background** — matte black fabric, black felt, or black grow bags produce the highest contrast against green plant tissue and are strongly preferred.
- Ensure the background material is **non-reflective**. Shiny or glossy surfaces can produce specular highlights that the algorithm may misclassify as plant tissue.
- The background should be **consistent across all images** in an experiment. Changes in background material between imaging sessions will introduce inconsistencies in the size measurements over time.

Avoid

- **Brown or tan substrates** such as exposed soil, coco coir, or peat visible in frame. These materials contain warm color tones that partially overlap with plant hue values and will be incorrectly included in the mask.
- **Green or blue-green containers**, trays, or pots. Any container color in the green hue range will be detected as plant material.
- **Patterned or textured backdrops** such as woven fabric with a visible weave pattern. Fine texture can produce noise at the segmentation boundary.

3.2 Substrate and Container Visibility

The growing substrate — soil, perlite, coco coir, rockwool, or similar media — should not be visible in the image frame. Even small exposed areas of substrate around the base of the plant can produce false positives in the mask, inflating the reported plant size.

Recommended

- **Cover exposed substrate** around the plant base with black fabric, black landscape cloth, or a fitted cover cut to the container opening. This eliminates the substrate from the camera's field of view entirely.
- Position the camera so that the **container rim is not visible**, or ensure it is a dark neutral color that does not fall within the green hue detection range.

- If using net pots or open containers, **wrap or cover the pot exterior** with black material before imaging.

Avoid

- Leaving **bare substrate exposed** around the plant stem. This is the most common source of segmentation error in practice.
- Using **white or light-colored grow media** such as perlite or hydroton without covering, as these reflect light and create bright regions the algorithm was not designed to exclude.
- **Wet substrate surfaces**, which can produce dark reflective patches that interfere with the background detection step.

3.3 Lighting

The segmentation algorithm performs best when the plant is **evenly and adequately illuminated**. The example images shown in this manual were captured under relatively low ambient light, which required the algorithm to perform additional brightness correction before segmentation. While the software handles underexposed images, consistent and diffuse lighting reduces processing error and improves measurement repeatability.

Recommended

- Use **diffuse, even illumination** across the plant canopy. Softbox lighting or LED grow panels with a wide beam angle reduce hard shadows between leaves.
- Ensure lighting is **consistent between imaging sessions**. Changes in light intensity or angle over the course of an experiment will cause apparent changes in plant size that are artefactual rather than biological.
- If imaging under natural or greenhouse light, **image at the same time of day** for every session to minimize variation in sun angle and intensity.

Avoid

- **Hard directional lighting** that casts strong leaf shadows onto the background or onto other leaves. Shadow regions on the plant can fall below the detection threshold and appear as holes in the mask.

- **Backlighting** or light sources positioned behind the plant relative to the camera, which reduces apparent leaf contrast.
- **Mixed color temperature light sources** in the same frame, for example a cool white grow light combined with warm incandescent ambient light, which shifts hue values unpredictably.

3.4 Plant Suitability

The segmentation algorithm detects plant pixels by identifying green hues within a defined range. Plants that fall outside this range at any point during the experiment will not be accurately segmented.

Compatible Plant Types

- Plants with **predominantly green foliage** throughout the experimental period.
- Species with **flat or broad leaves** viewed from above produce the cleanest masks, as leaf area is maximally exposed to the camera.
- Plants at **vegetative growth stages** where green leaf tissue dominates the canopy.

Incompatible or Problematic Cases

- **Severely stressed or senescing plants** where significant yellowing, browning, or necrosis has occurred. Yellow tissue falls outside the green hue detection range and will be excluded from the mask, underreporting true plant size.
- **Purple or anthocyanin-rich varieties** such as red lettuce or purple basil. These cultivars contain little green pigment and will produce very small or empty masks.
- **Seedlings at germination stage** where only the hypocotyl and cotyledons are visible. These structures are very small and may fall below the minimum area threshold used to filter noise. The fill threshold parameter described in the Settings subsection of the Hardware and Software chapter (Section [1.2.4](#)) may need to be reduced for early-stage imaging.

- Plants where **flowers, fruit, or senescent tissue** constitute a large fraction of the visible canopy, as these structures are typically non-green and will be excluded from the size measurement.

3.5 Pre-Experiment Setup Checklist

Before beginning an imaging experiment, verify each of the following conditions:

- ☐ Background is dark, matte, and uniform across all imaging positions.
- ☐ All exposed substrate around the plant base is covered with black material.
- ☐ Container exteriors are dark or covered and do not fall within the green hue range.
- ☐ Lighting is diffuse, even, and will remain consistent across all imaging sessions.
- ☐ Plant species is majority-green and at a vegetative growth stage.
- ☐ Camera position is fixed and will not shift between imaging sessions.
- ☐ A test image has been processed through the Analysis function and the resulting mask reviewed before committing to a full experiment run.

Chapter 4

Troubleshooting

Problem	Resolution
No Excel files found during preview	Verify that the correct root folder has been selected and that the Analysis step has been completed for at least one plant.
Error processing a specific file	Check that the file is not open in another application. Ensure the file was generated correctly by the Analysis function and has not been manually modified.
Merged.xlsx is missing charts	Ensure that plant groups have been assigned before clicking Merge Sheets. If groups are unassigned, all plants default to the Control group color.
Matplotlib charts are not generated	Verify that the Generate Matplotlib Plots checkbox is enabled and that at least one plant is not excluded in the Exclude Plants window.
Software becomes unresponsive during processing	This may indicate insufficient system memory. Verify that available RAM meets the minimum requirement of 12 GB and close other applications before restarting the software.
Debug window shows output from the wrong tab	Switch to the desired tab before running any operation. The debug output window updates to the active tab automatically.

Chapter 5

Future Considerations

The SMART Platform Imaging Software has been designed to provide an extensible platform for automated plant imaging and data analysis. Future development may include improvements in both the hardware and software to increase functionality, usability, and processing efficiency.

5.1 Hardware

Potential hardware enhancements include:

- Installation of limit switches at the maximum travel extents to provide additional protection against out-of-bounds movement.
- Integration of improved lighting systems to provide more consistent image quality across experiments.
- Support for higher-resolution or industrial machine vision cameras.
- Light absorbing background or increased layering of landscape fabric to decrease luminance in images.

5.2 Software

Potential software enhancements include:

- Automatic software update functionality.
- Support for additional image analysis algorithms and plant phenotyping metrics.

- GPU acceleration to reduce processing time for large datasets.
- Automatic detection of hardware devices, eliminating manual configuration of serial ports and camera indices.
- Expanded data visualization and statistical analysis tools.
- Additional export formats for processed data and reports.
- Cross-platform compatibility with Linux and macOS.
- Support for various varieties of various color ranges — not just primarily green varieties.

The modular design of both the SMART platform hardware and the SMART Platform Imaging Software allows new capabilities to be incorporated with minimal changes to the overall workflow.